



spectrophotometers and test kits

Spectroquant® Prove for accurate and convenient water, beverage & environmental analysis



The life science busines of Merck operates as MilliporeSigma in the U.S. and Canada.

Supelco®
Analytical Products

spectroquant® prove

As your trusted analytical partner, we built the perfect tool for water analysis, one that unites the simplicity you want with the confidence you need, and the durability you expect. Spectroquant® Prove delivers all this and more. Its intuitive design and our comprehensive portfolio of Spectroquant® test kits make analysis smoother than ever and support you throughout your workflow.

Long-life Lamp

Extended durability and economy

The new reference beam technology in our Spectroquant® Prove instruments minimizes the effort and costs of replacing lamps. Hence, the halogen lamp in the Prove 100 lasts significantly longer, while the special xenon lamp in the Prove 300 and 600 offers a lifetime of several years combined with the flexibility of UV/Vis analysis - even under intense use.

Smart Screen

Simple, seamless menu navigation

Have complete control at your fingertips with our interactive Smart Screen and intuitive icons.

- Method and result lists for easy documentation and data access
- Touchscreen is resistant to most chemicals
- System menu in 28 built-in languages

Custom Measurements

Easy adjustments for every need

Several modifications can be made easily using the settings icon, for example:

- Simultaneous display of absorption and concentration values
- Turbidity correction
- Add the dilution factor for each sample

Individualized Measuring Ranges

Set your own individual ranges:

- Monitor results within specified limits
- Visualize expected ranges
- Simplify your trend analysis



Data Exchange

Fast and flexible connectivity

You have multiple ways to easily utilize your data:

- USB and Ethernet port to connect your printer or memory device
- "Prove Connect to LIMS" allows you to push your results directly to your LIMS or ELN (Cat.No. Y1100860001)



Ready-to-use Test Chemistry

Maximum flexibility

Every Prove model accepts 16 mm round cells as well as 10, 20 and 50 mm rectangular cells (cuvettes). For even greater sensitivity, the Prove 600 is compatible with 100 mm pathlength cells.



Minimal effort

The Prove cell test port offers direct insertion of our readyto-use round cell tests without opening any lids. The Live ID barcode identifies the correct method instantaneously.

Our innovative AutoSelector technology means that your reagent testing will be precise and consistent. It recognizes the cell size, detects the analysis method and calculates the results based on the cell size.



Easy Cleaning

Durable, removable components

The intelligent design makes it easy to keep your Spectroquant® Prove in top form – even if something is spilled into the sample compartment. Simply remove the cell holder, rinse it with water, wipe it clean, and it's ready for your next analysis.

Compact

Dimensions: W 42 cm, D 28 cm, H 24 cm

Accuracy and precision made simple

Live ID

A new dimension in traceability


Live ID is an identification system for all our cell and reagent tests. Unlike other barcodes, it transfers "live" data—thanks to automatic calibration updates, meaning that your results will always be accurate and fully traceable. Live ID also saves time; just insert the test, and your Prove instrument reads the code, detects the method, and selects the right settings. Test-specific lot and expiry data are saved with each result for accurate records and streamlined audits.



- Method recognition
- Lot number
- Expirv date
- Calibration update

Reagent Tests with AutoSelector

These economical kits contain highly stable reagent mixtures, and allow you to choose the cell format that suits your measuring range. An AutoSelector with a Live ID barcode is provided with each reagent test kit for accurate measurements.

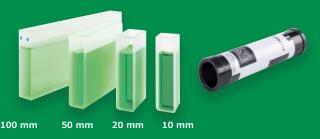



Cell Tests

For maximum convenience, our cell tests include all the reagents you need for your analysis, so they are ready to use straight out of the box.

Spectroquant® Test Kits

Not just more choices. More confidence.

Spectroquant® Prove is preprogrammed for over 180 Spectroquant® test kits and methods, giving you swift, reliable results for the broadest range of applications.





Data on Demand

Customizable data for easy documentation

- Enjoy maximum flexibility in minimum steps
- Access all your results from a clear, concise list
- Select your desired filter option (by date, keyword, method, user, etc.)
- Export data in your preferred form from brief overview to complete protocol



AQA Prime

A new level of analytical standard

AQA Prime makes analytical quality assurance exceptionally reliable - and extremely easy.

- Clear overview of AQA status, information, updates, and options with one touch
- · Select individual settings and measurements for each AQA
- Complete verification and documentation for all stages of quality control:

AQA1 > Instrument check

AQA2 > System check – now even easier because Spectroquant® Prove is preprogrammed for all Spectroquant® standards













Matrix Check

Easily detect interferences

The Matrix Check makes it even easier to spot foreign substances in your samples. After your measurement, simply press the Settings icon, and select Matrix Check from the list. You will be clearly guided through the entire process.

Preprogrammed Applications

Free, additional methods to simplify your work

To further facilitate your analysis, Spectroquant® Prove offers free, preprogrammed application procedures, such as bromate in drinking water, ICUMSA and brewery methods.

Robust Surfaces

Built to endure lab chemicals

Spills and splashes happen no matter how careful you are. That's why the materials and design of Spectroquant® Prove are engineered to last and be resistant to most chemicals used in labs.

spectroquant® prove

100 | 300 | 600

Spectroquant® Prove 100

Cat. No. 1.73016.0001



For routine applications

The Prove 100 is the best choice for those who primarily use our broad range of Spectroquant® test kits, or only perform Vis measurements. High quality and great value for your daily analyses.

Spectroquant® Prove 300

Cat. No. 1.73017.0001



For sensitive measurements

Thanks to its long-lasting xenon lamp, the Prove 300 is ideal for more intensive use. What's more, it is capable of both UV and Vis measurements, so you have greater flexibility for more intricate analyses.

Spectroquant® Prove 600

Cat. No. 1.73018.0001



For complex analyses

Designed with high-end UV/Vis optics, and cuvettes of up to 100 mm, the Prove 600 packs great power into a compact size. Excellent resolution and sensitivity for complex kinetics or spectral measurements, or for use with test kits.

Specification	ns	Prove 100	Prove 300	Prove 600
Measuring technology	Spectrophotometer with reference beam technology			
Wavelength range	Vis (320 - 1,100 nm)			
	UV (190 - 320 nm)	•		
Lamp type	Tungsten halogen lamp			
	Xenon flash lamp			
Ambient light protection	Measurement with open shaft possible due to proprietary solution (patent pending)			
Spectral bandwidth	4 nm			
	1.8 nm			
Smart Screen Display	Resistive touchscreen			
	P-cap glass touchscreen			
Live ID System	2-D Barcode recognition for cell tests and reagent tests			
	Barcode contains lot, expiry, and calibration data; data stored with each measurement			
Cell size	16 mm round cells, 10, 20 and 50 mm rectangular cells with automatic recognition			
	100 mm rectangular cells with automatic recognition			

Specification	s	Prove 100	Prove 300	Prove 600
Cell holder	Removable for easy cleaning			
Methods	Programmed methods of all Spectroquant® cell and reagent tests, 99 user-defined methods, 20 profiles for kinetics and absorbance scans			
Applications	Free preprogrammed applications: Bromate, Brewery packages (MEBAK/ EBC methods), Sugar (ICUMSA), Oil (DOBI, olive oil)			
AQA Prime	Individual settings for all methods in AQA 1 (instrument check) and AQA 2 mode (system check) and pipette check			
Sample matrix check	Instrument-supported matrix check for each method			
Software updates	Free method updates on the MilliporeSigma website			
Languages	Navigation in 28 built-in languages			
Communication interfaces	USB: 2 x USB-A (for printer, USB memory devices, keyboard or barcode reader), 1 x USB-B; Ethernet: LAN connection			

Application fields



Wastewater



The Prove 100 is the perfect choice for routine Vis measurements. Add your wastewater sample to the cuvette or cell, then insert and read your results. The Spectroquant® portfolio offers a broad selection of COD and other test kits in multiple measuring ranges.

Drinking water & beverages



The analysis of drinking water and beverages typically involves reagent tests since they offer lower detection limits for parameters such as manganese and sulfate.

The Prove 300 is ideal here, as it allows UV/Vis analyses, and is preprogrammed with free applications, like bromate and brewery tests.



Process water



Even low levels of impurities in process water can lead to damage, downtime, and expensive repairs. To help you avoid these, we have developed the most sensitive silicate and chloride tests with the lowest measuring ranges in the market. For even greater sensitivity, use the 100 mm cell in the Prove 600.

Application notes

Our photometric application notes demonstrate a broad range of methods to support your production and analysis.

Whether you need to determine COD in wastewater, the level of nitrate in soils or the concentration of bromate in your drinking water - over 200 application notes support your sample preparation and analysis.

Compliance, always ensured.

At Merck, we understand the importance of reliable instruments, and the complexities of documentation. That's why we offer a comprehensive service plan that takes care of both.

All Spectroquant® service plans offer the following features:

- Performance check with reference materials certificate included
- Factory-recommended maintenance to prevent breakdowns
- Telephone support hotline connects you to our technical experts
- Free software upgrades to keep your instrument up-to-date
- Reagent delivery program for convenient, punctual supply of our high-quality test kits

Supelco® Analytical Products



© 2019 Merck KGaA, Darmstadt, Germany and/or its affiliates. All Rights Reserved. Merck, the vibrant M, Supelco and Spectroquant are trademarks of Merck KGaA, Darmstadt, Germany or its affiliates. All other trademarks are the property of their respective owners. Detailed information on trademarks is available via publicly accessible resources.

Lit. No. MK_BR1091EN Ver. 1.0 2018 - 18838 08/2019

www.scientificlabs.co.uk

